


# ALL IN ONE MEDIA PREPARATION


**LABGIG**



## SIMPLIFY

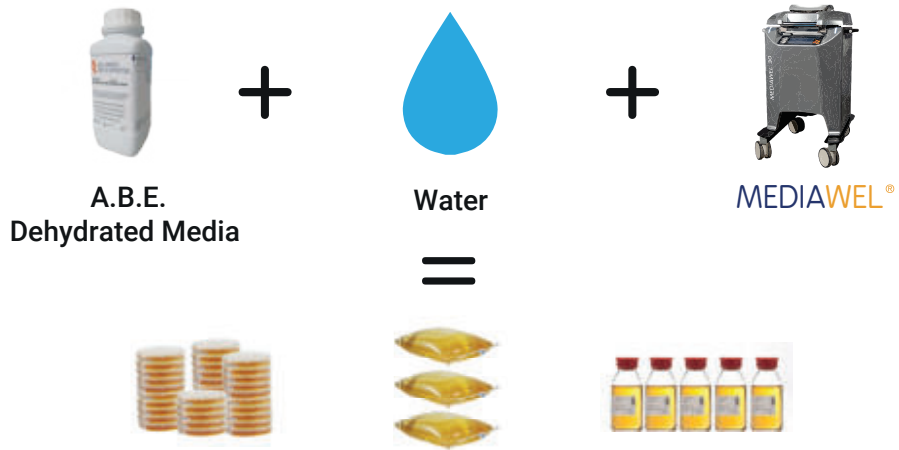


## HOMOGENEITY

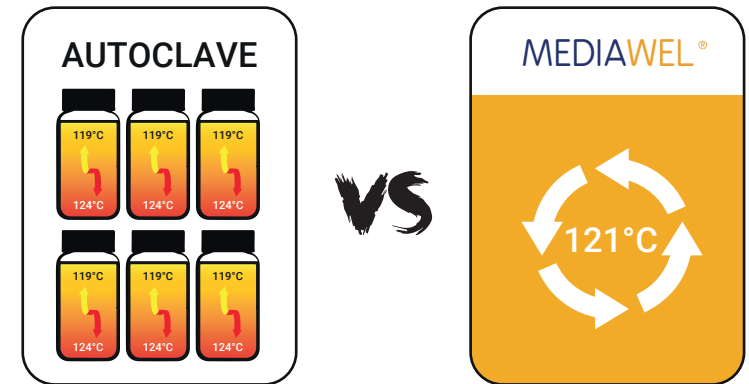


## ALL IN ONE MEDIA PREPARATION

## TIME SAVING

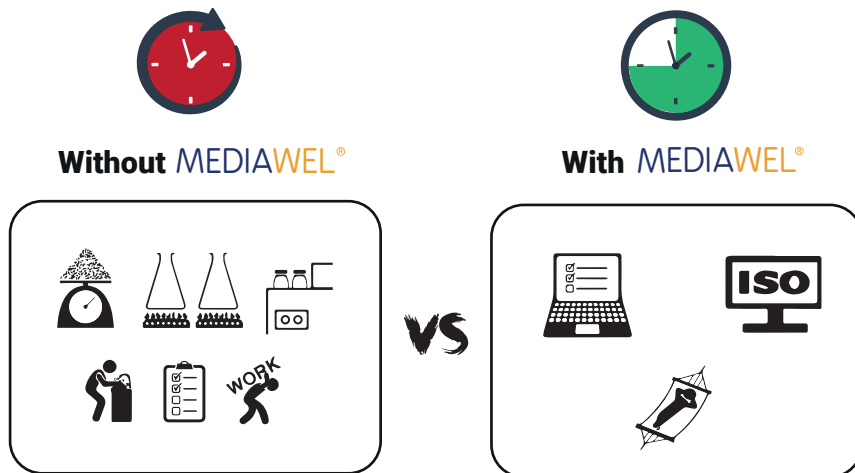

## FAST COOLING SYSTEM

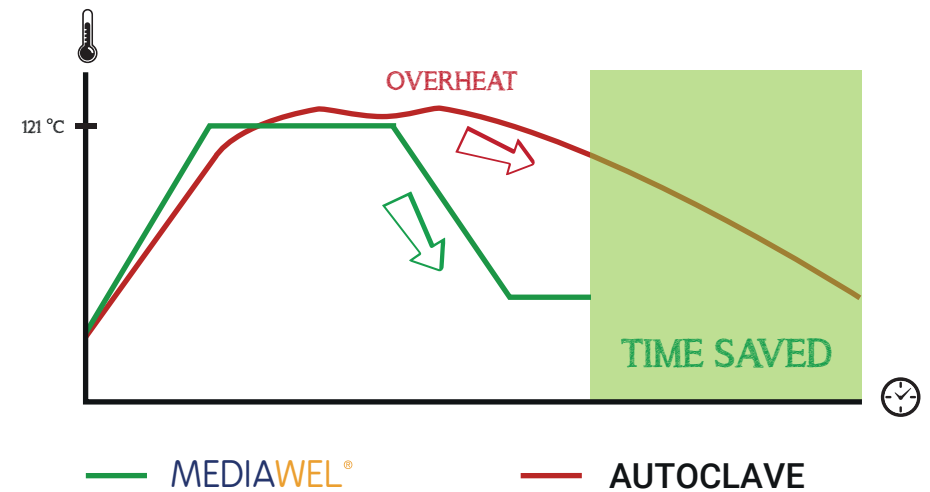


# MEDIA PREPARATOR

Prepare up to 50L of media

# MEDIAWEL 50<sup>®</sup>



**Safety lid**  
Locked above 80°C

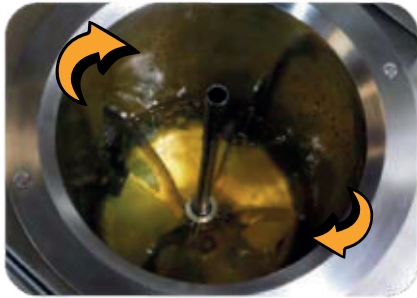

**Large stirrer**  
From liquid to viscous media



**e-TRACE<sup>®</sup> Technology**



**Easy to move**  
Large wheels with brakes



**Multicolored visual and audible alarm**

**7" LCD touchscreen**  
Gloves friendly

**Antishock body**  
ABS PC Plastic



**Optimized compacity**  
Same footprint as a 30L media preparator :  
MEDIAWEL 30<sup>®</sup>

## PERFORMANCES

From 5 to 50 liters

60 - 90 min

25-125 °C / ± 0.5 °C

**Delivered with MEDIAWEL 50** : USB stick with user Manual, control certificate, large magnetic stirrer (MEDW2008), sampling tube (MEDW2009), draining plug (MEDW2028), sampling nozzle seals x2 (MEDW4011), sampling nozzles x2 (MEDW2003 : Ø 12mm + MEDW2004: Ø 6.4mm), inlet and outlet water tubing (4 m) + 4 tubing clamps (MEDW4012), decompression filters x2 (MEDW3001), Bactystopper - calibrated sterile carded cotton x50 (MEDW3003), sliding washer stirrer axis x2 (MEDW4003), 3 Months free access to the Labpage "Real time data monitoring LICENCE"

## TECHNICAL DATA SHEET


Capacity	From 5 to 50 liters
Sterilization temperature	From 25°C to 125°C
Sterilization time	From 1 to 99 min
Distribution temperature	From 25°C to 80°C
Temperature accuracy	+/- 0.5°C after temperature stabilization
Remaining volume	<130ml
Adding cap	✓
Traceability	Integrated eTRACE® module, USB export of PDF/CSV files
WiFi/Ethernet traceability	Labpage® dashboard. Real time data monitoring, data storage on SD card (up to 10 000 cycles)
Temperature probe	3 probes (1 standard, 1 in option and 1 for calibration)
Safety pressure release	1.7 bar safety valve
Dispensing port	Adapters for tubes of 6.4 mm / 12 mm
Continuous homogenization	✓
Touchscreen	Tactile, LCD 7"
Alarm	Multicolored visual and audible alarm
Cleaning	Bottom drainage / cleaning evacuation
Mobility	Stable and easy to move : 4 large medical wheels with brakes
Robustness	Antishock body (ABS PC PLASTIC)
Vessel	Stainless steel AISI 316L
Programs	50
Weight	178kg
Dimensions (L x l x H)	735 x 585 x 1073
Warranty	1 year
Power supply	400v +/- 10% - 50/60Hz - 28A (Max by phase) 17000 W
Country of origin	France

# MEDIAWEL 50®



# TRACEABILITY WITH YOUR MEDIAWEL®

Each cycle is exported (PDF/CSV) to your USB key :

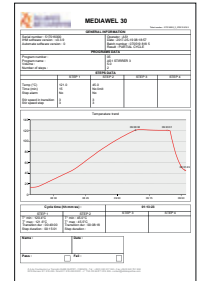

**REAL TIME DATA MONITORING** Labpage® Real time data monitoring through WIFI/Ethernet, connected to your PC / Smartphone / Tablet :

**Dashboard for the live monitoring :** program settings, temperature settings and current temperature, alarms and notifications ...



Cycles are automatically recorded and stored. Recall them one by one or sort them out by batch number, date, user ...



\*SD card included - up to 10 000 cycles data storage